
© 2014 Wahl Clipper Corporation  
Part No. 93342-4801 Printed in China

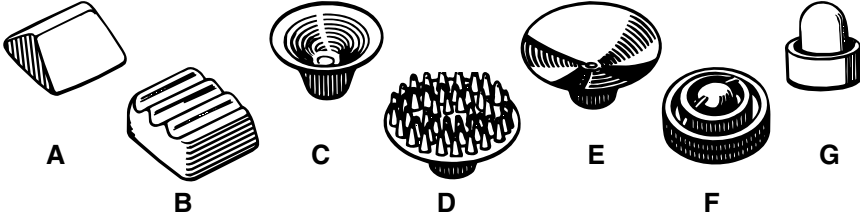WAHL CLIPPER CORPORATION  
World Headquarters  
P.O. Box 578, Sterling, IL 61081

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Wahl (UK) Ltd.  
Hemry Trade Park  
Sea Street, Herne Bay, Kent  
CT6 8JZ, England

Wahl Clipper (PTY) Ltd.  
Suite L9 & L10  
Hingham Field Office Park  
79 Boeing Road East  
Bedfordview, Gauteng 2007, South Africa

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Wahl India  
PO Box 456  
Acacia Ridge QLD 4110  
Australia

Wahl India Grooming Products Pvt Ltd  
201, 2nd Floor, X Cube, Opp. Fun Republic  
New Link Road, Andheri (West)  
Mumbai – 400 053 Maharashtra, India

Please visit [www.wahlglobal.com](http://www.wahlglobal.com)

## en Operating instructions

### IMPORTANT SAFEGUARDS

For your own protection from injuries and electric shocks, the following information must be observed when using electrical equipment:

- Read all the instructions before using the appliance!
- Never touch electrical equipment that has fallen into water. Disconnect from the mains immediately.
- Never use electrical equipment in the bath or shower.
- Always position or store electrical equipment where it cannot fall into water (e.g. wash basin). Avoid electrical equipment being exposed to water and other liquids.
- Always disconnect the appliance from the power supply immediately after use except when charging.
- Disconnect the appliance from the power supply before cleaning.
- This appliance is not suitable for use by people (including children) with physical, sensory or mental disabilities or who do not have the relevant experience or knowledge of using such appliance unless under the supervision of or after having received instruction in the use of the appliance by an individual who is responsible for their safety. Children should be supervised to ensure they do not play with the appliance.
- Only use the appliance for the proper purpose as described in the operating instructions. Only use accessories recommended by the manufacturer.
- Never use the appliance if the mains cable or plug-in transformer is damaged. In addition, do not use the appliance if it is not working correctly, if it is damaged or if it has fallen into water. In these cases, send the appliance to our Service Centre to be inspected and repaired. Only specialists properly trained in the handling of electro technical devices may repair electrical equipment.
- Do not hold the appliance by the mains cable to carry it and do not use the cable as a handle.
- Keep the appliance away from hot surfaces.
- Make sure that the mains cable is not twisted or kinked during charging or stornig.
- Never insert or allow objects to fall into the openings of the appliance.
- Only use and store the appliance in dry rooms.
- Never use the appliance where aerosol (spray) products are being used or where oxygen is released.
- The battery must be removed from the appliance before it is scrapped. The appliance must be disconnected from the supply mains when removing the battery. The battery is to be disposed of safely.
- Cleaning and user maintenance shall not be made by children without supervision.
- The appliance is only to be used with the power supply unit provided with the appliance.
- Keep the appliance dry.
- The appliance is safety insulated and radio screened. It meets the requirements of EU Directive 2004/108/EC regarding electromagnetic compatibility and EC Low Voltage Directive 2006/95/EC.

Do not use massager without attachments.

- To reduce the risk of burns, fire, electric shock or injury to persons:
  - This appliance should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.
  - Do not operate this appliance under a blanket or pillow. Excessive heating can occur and cause fire, electric shock or injury to persons.
  - In case of pregnancy, diabetes or illness, consult your physician before using this massager.
- Do not use massager on:
  - any sensitive skin, swollen or inflamed areas, in the presence of poor circulation, where skin eruptions are present, or in the presence of unexplained calf or abdominal pain.
  - Varicose veins.
  - The front of the neck.
  - Genital areas.
  - Sleeping or unconscious persons.
  - When using a massage attachment, be sure to seat the attachment fully onto the stern and keep it on while massaging.
  - Do not use this massager for more than 20 minutes at one time.
  - Always attach plug to appliance first, then to outlet.
  - Plug charger directly into receptacle--do not use an extension cord.
  - To disconnect, turn unit to the OFF position, then remove plug from outlet.

RECHARGING

- This massager is designed for household use. It should be allowed to charge 19 hours before first use.
- Connect the transformer to the massager. BE SURE MASSAGER IS IN THE "OFF - 0" POSITION WHEN CHARGING.
- Plug the recharging transformer into any convenient electrical outlet at the correct voltage.
- Ni-MH battery life may be extended with the following procedure: Once every month, let the massager run until the motor stops. Move the switch to the "OFF - 0" position and recharge fully overnight.
- The charging indicator light will glow when your massager is in the "OFF" position and connected to its charger. The massager cannot be overcharged, so you may charge it as often as desired. We recommend that the battery be "worked" or discharged by massaging without the cord and by not recharging after every use. In certain areas that experience low line voltages, the red charging indicator light on your massager may become dim as the battery nears full charge. This dimming does not affect the performance of your massager but is simply a result of variances in power line voltages.

### TROUBLESHOOTING

- Your massager does not seem to take a charge, check the following before returning it for repair:
  - Check outlet current by plugging in another appliance.
  - Make sure the outlet simply is not connected to a light switch which turns the power at the outlet off when light is off.
  - Be sure the massager is properly connected to the recharging transformer, all connections are secure and then switch on the massager is turned "OFF - 0".

### OPERATING INSTRUCTIONS

- Place the desired massage attachment on the stern of the massager. The attachments easily press on and pull off.
- Apply the massager to the targeted area, using continual movement. A circular motion is often most effective. (See MASSAGE TECHNIQUES.)

### MASSAGE ATTACHMENTS

- KNUCKLE JOINT ATTACHMENT (A)**
  - For concentrated action between the fingers or toes. Great for knee and elbow joints, too.
- MUSCLE KNEADING ATTACHMENT (B)**
  - Delivers ripple-effect action for large muscle areas.
- GENERAL BODY ATTACHMENT (C)**
  - For general, all-body massage.
- SCALP STIMULATOR (D)**
  - Delivers effective scalp massage, loosens dandruff, and is soothing and relaxing.
- FacIAL VITALIZER (E)**
  - Provides beneficial, revitalizing facial massage. (Apply cold cream or lotion to face prior to use. Keep massager away from any eyes.)
- DEEP MUSCLE STIMULATOR (F)**
  - Targets deep into tissue. For general, all-body massage.
- SPOT APPLICATOR (G)**
  - Great for pinpointing concentrated areas.
- CHANGING ATTACHMENTS**
  - To remove an attachment head, use your thumb and two fingers, and pull away from the massager stem. To insert an attachment head, align the head with the massager stem, and press firmly onto stem.
- CLEANING/STORAGE**
  - Always unplug massager before cleaning.
  - To clean, simply wipe unit with a dry, soft cloth. Do not immerse in water or use any harsh cleaners or solvents. When storing unit, make sure cord is not twisted or kinked.

### MASSAGE TECHNIQUES

- SHOULDERS AND NECK** — Massage each shoulder toward back of neck to point just behind and below ear. Then back to shoulder. Repeat several times.
- ARMS** — Start with a gentle massage of the palms and fingers. Next, using a circular motion, move upward to forearms. Then turn up to biceps and shoulders.
- LEGS** — If desired, massage bottoms of feet. Then, using a circular motion, massage each leg from calf on up and over thigh.
- CHEST** — Massage from center of chest outward, and up toward each shoulder.
- BACK** — Massage from top of spine to bottom, and from center outward.
- ABDOMEN** — Relax abdominal muscles by massaging from center outward, using a gentle, circular motion.
- JOINTS** — Relieve tired, aching knees and elbows with smooth, circular strokes.
- FACE** — Gently massage from the nose outward, taking care to avoid contact with the eye.
- SCALP** — Start at the back of the neck and massage up to the top of the head. Using this up and down motion, massage across the scalp.

### Disposal in EU countries

Do not dispose of the equipment with household waste. As part of the EU Directive governing the disposal of electric and electronic appliances, the appliance is accepted free of charge by local waste collection points or recycling centres. Correct disposal will ensure environment protection and prevent potential harmful impact on people and the environment.

### Disposal in non-EU countries

Please dispose of the equipment at the end of its life in an environmentally friendly manner.

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